


DIAGRAMMATIC RECAPITULATION OF APPLICATION NOTE NO. 14

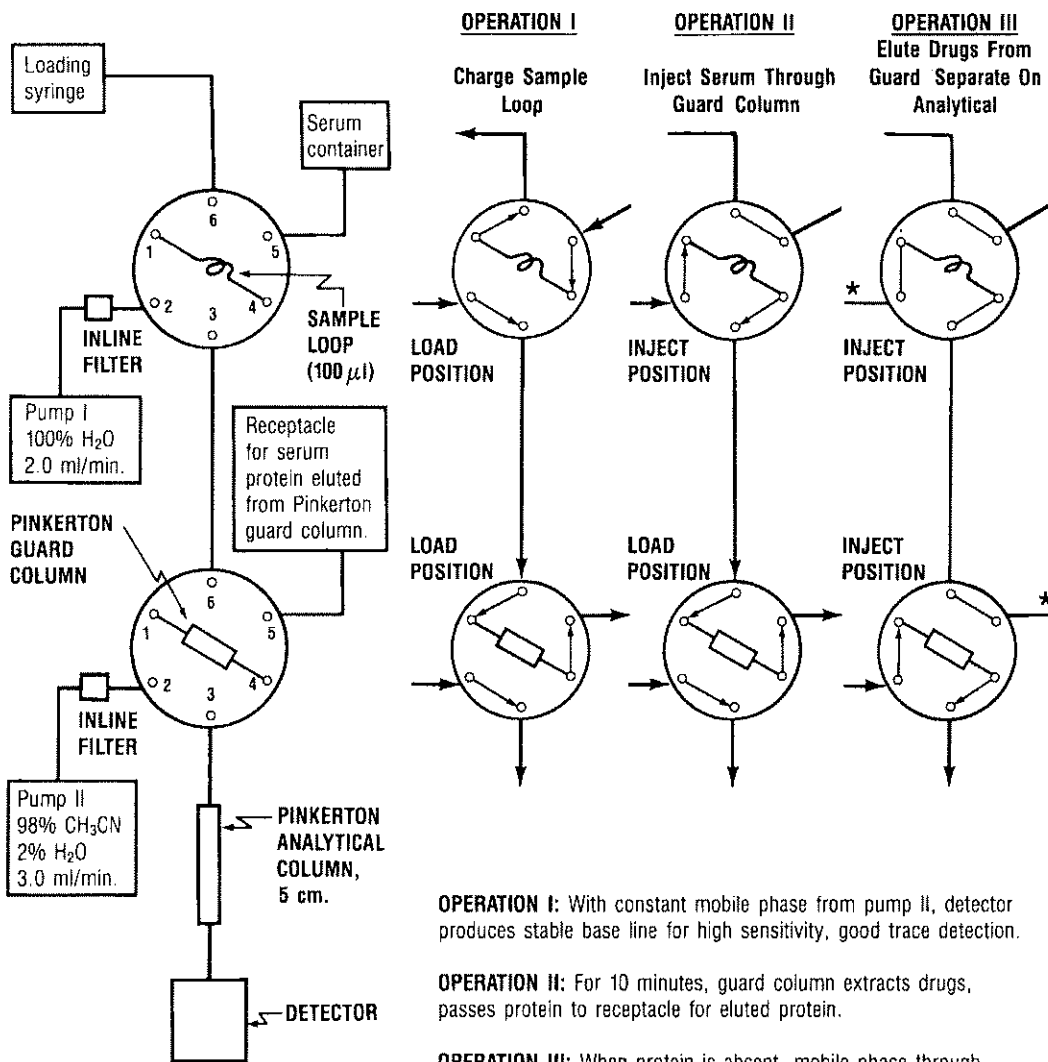

VALVE: REGIS CODE 701000
 SAMPLE LOOP: REGIS CODE 701010
 INLINE FILTER: REGIS CODE 731407

OPERATION I: With constant mobile phase from pump II, detector produces stable base line for high sensitivity, good trace detection.

OPERATION II: For 10 minutes, guard column extracts drugs, passes protein to receptacle for eluted protein.

OPERATION III: When protein is absent, mobile phase through Pinkerton analytical column may be 98% CH₃CN, as it now is.

*Pump I stopped.

T.J. Szczerba, Regis
 Chemical Company, 8210
 Austin Avenue, Morton
 Grove, IL 60053